

A
Monteringshuller i gulvet - $\varnothing$ 10 mm.
Trous de montage - $\varnothing$ 10 mm.
Montierungsslöcher - $\varnothing$ 10 mm.
Mounting holes - $\varnothing$ 10 mm.

B
Koldt vand - R 3/8"
Amenée d'eau - R 3/8"
Wasserzufluß - R 3/8"
Water inflow - R 3/8"

C
Trykluft - R 3/8"
Air comprimé - R 3/8"
Preßluft - R 3/8"
Air inflow - R 3/8"

D
Net 230 volt - 3x1.5
Secteur 230 volt - 3x1.5
Netz 230 volt - 3x1.5
Power supply 230 volt - 3x1.5

E
Klokkekald - 2x1.5
Ligne de sonnerie - 2x1.5
Klingelleitung - 8x0.5
Bell lead - 2x1.5





F
Styreledning til sug - 2x1.5
FConducteur pour pompe aspirante - 2x1.5
Schaltleitung für Saugmaschine - 5x1.5
Control lead for suction machine - 2x1.5

G
Røntgenapparat - 6x1.5 afhængig af fabrikat.
Radio scopique - 6x1.5 dépendent du produit fabriqué
Röntgenapparat- 6x1.5 abhängig von der Fabrikat.
X-ray equipment - 6x1.5 depending on the manufacture.

H
Sug - $\varnothing$ 32 mm
Conduit d'aspiration - $\varnothing$ 32 mm
Saugleitung - DN40
Suction pipe - $\varnothing$ 32 mm

J
Vandafløb - $\varnothing$ 32 mm
Evacuation d'eau sans filetage - $\varnothing$ 32 mm
Wasserabfluß ohne Gewinde - DN40
Water outflow without thread - $\varnothing$ 32 mm

	Instrument table height - A	Main arm height - B	Overall height - C
Low column	~810 - ~1120	1320	1743
Middle column	~915 - ~1225	1425	1849
High column	~1015 - ~1325	1525	1949

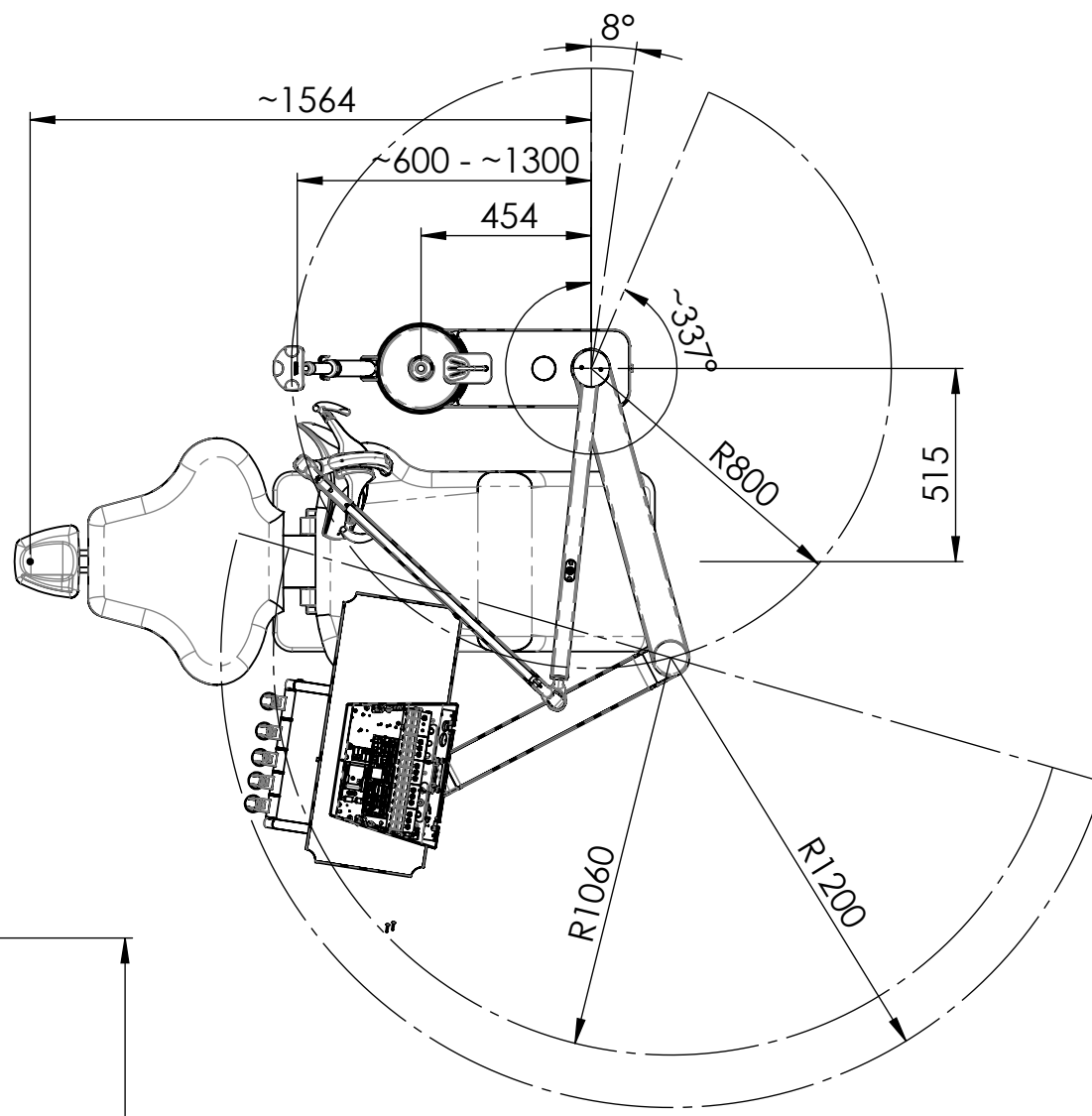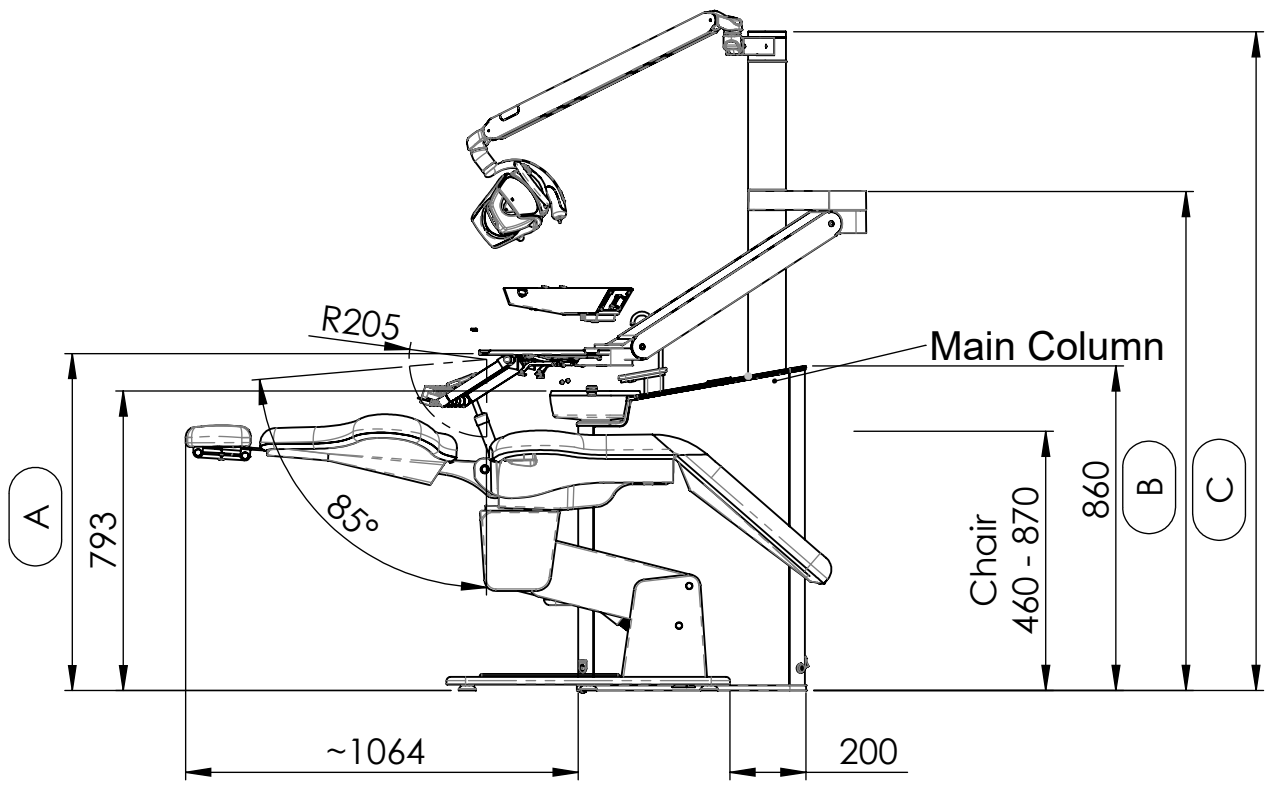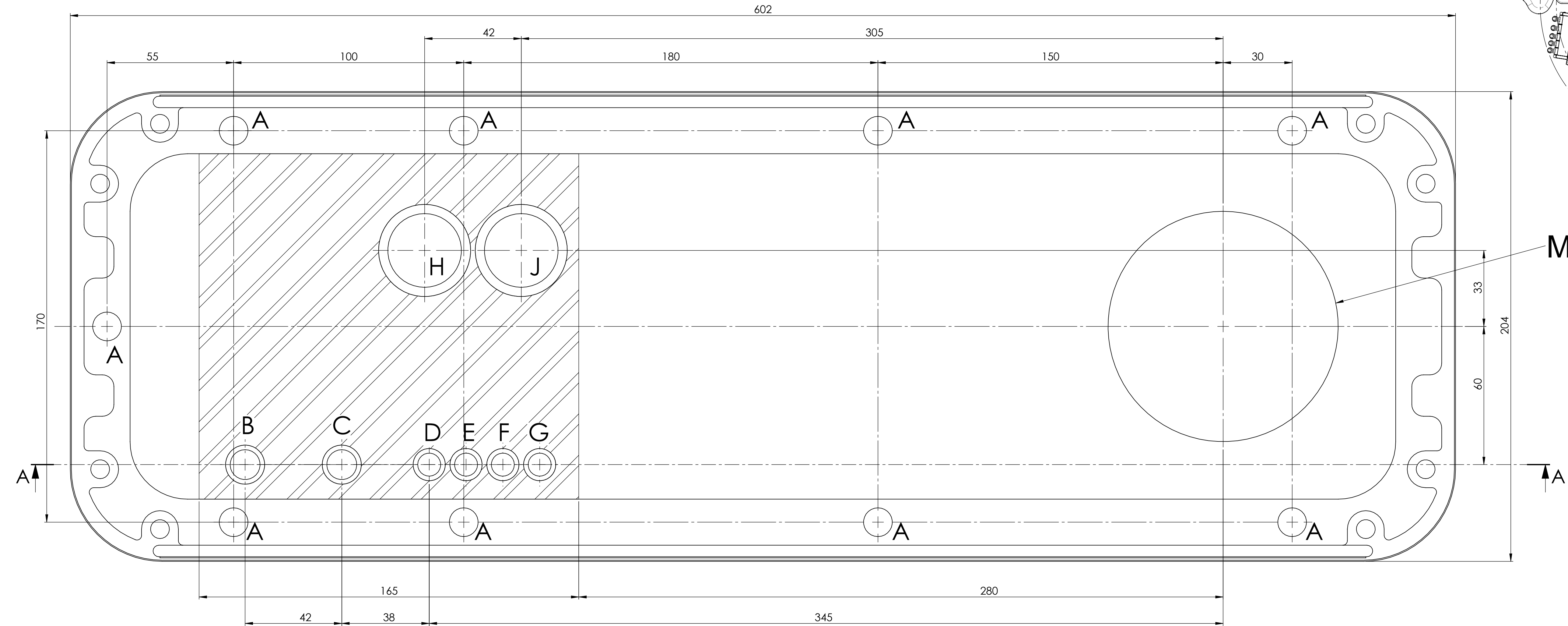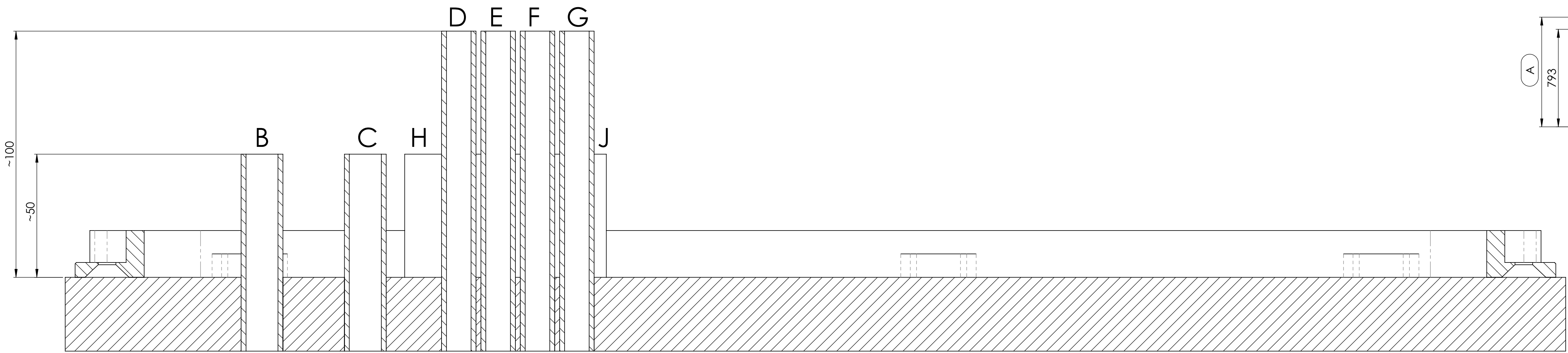

Main Column

HEKA DENTAL		Litauen Alle 4, DK-2630 Taastrup, Denmark Phone: +45 43 32 09 90, Fax: + 45 43 32 09 80 E-mail: mail@heka-dental.com	
Item name: Overview of Heka G+		Drawing no.: ar-0117-06	Rev.: 02
Scale: 1:1		Drw.by: pl 23/11-2023	

THE INFORMATION CONTAINED IN THIS DRAWING IS THE SOLE PROPERTY OF HEKA DENTAL. ANY REPRODUCTION IN PART OR AS A WHOLE WITHOUT THE WRITTEN PERMISSION OF HEKA DENTAL IS PROHIBITED.